


Mounting a Joystick

Here is a guide to the different mounting options and how to attach them.

Tabletop – Hook and Loop Fastener

One way the joystick can be mounted is on a tabletop or other surface, such as a lap tray, and fastened using a hook and loop fastener, such as Velcro.

Typically, best practice is to stick the side with the loops (the softer side) to the tabletop or other surface, and to stick the side with the hooks (rough side) to the joystick as shown below.

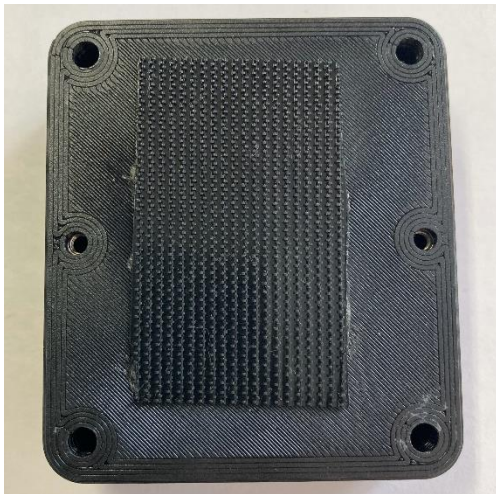

Tabletop – Non-Slip Pads

Another way the joystick can be used on the tabletop is by using a non-slip material, such as the no slip stoppers shown below, or a non-slip material such as Dycem.




Joystick Mounting Guide

¼-20 Threaded Camera Mount

Another way the joystick can be mounted, is by using a mount with a ¼-20 thread on the end, such as a commercially available camera mount. To mount the joystick on a ¼-20 thread, the Camera Mount Adapter is used.

Camera Mount Adapter:

When mounting on a mount with ¼-20 thread, the Camera Mount Adapter can be attached to a joystick using two M3 screws that are at least 10 mm long.

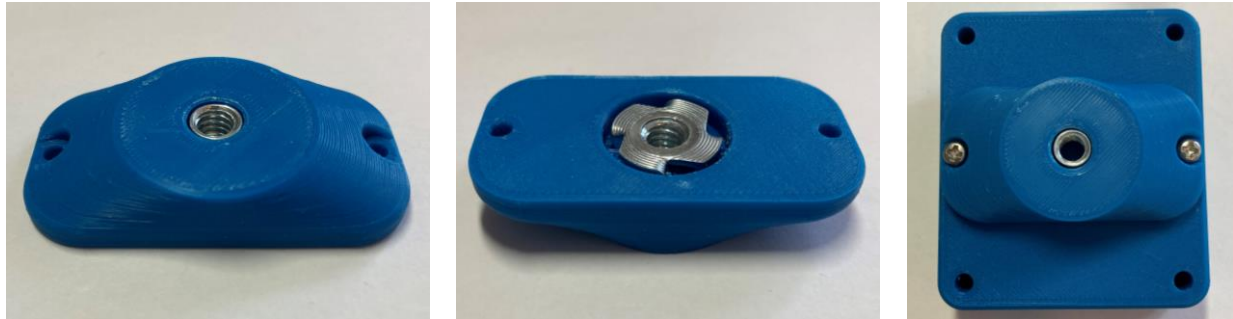

RAM Mount

Another way the joystick can be mounted is by using a commercially available RAM B compatible mount with a 1" ball and socket system. For this, the Joystick RAM Mount Adapter is used.

RAM Mount Adapter

When mounting on a RAM B compatible mount, the RAM Mount Adapter can be attached to a joystick using two M3 screws that are at least 10 mm long.

